


Holding Micropipets

From the First Name in Micropipets

- Humagen Holding pipets have a flat, polished opening with a large ratio of outer diameter to inner diameter for maximum control.
- Specially designed equipment ensures each pipet has a consistent outer diameter and flat opening.
- Large selection of standard sizes offered.

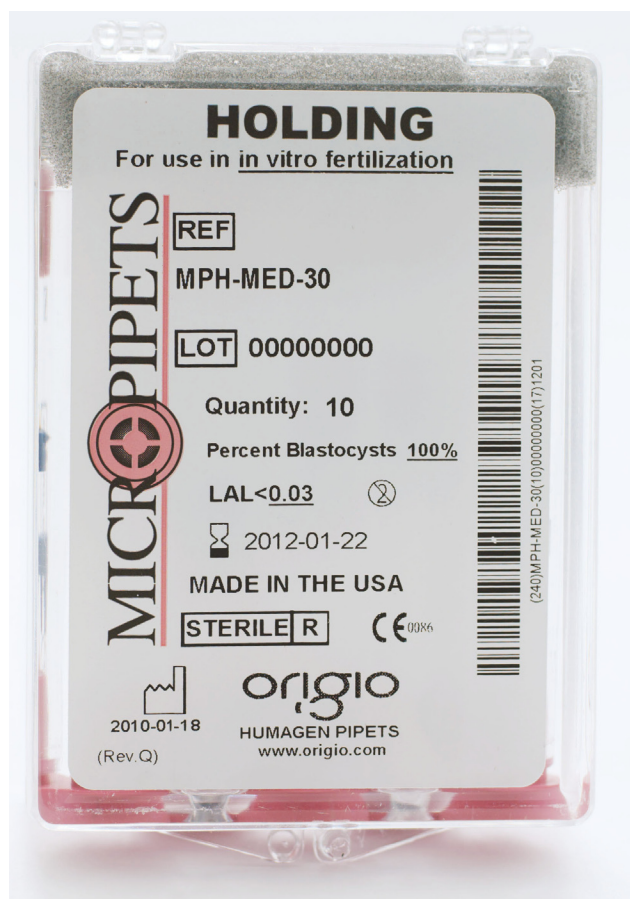


Each type of micropipet is color-coded for easy identification between pipet types.

Holding Micropipets

Holding micropipets are used to hold and immobilize the oocyte or blastocyst.

- Large and Extra large sizes are ideal for handling blastocysts.
- Custom sizes may be ordered by specifying the exact size (O.D. or I.D.) rather than the indicated ranges

Holding Micropipets (10 per box)

	Reference Number	I.D.	O.D.
	MPH-SM-(angle)*	15-20µm	65-95µm
	MPH-MED-(angle)*	15-20µm	95-120µm
	MPH-LG-(angle)**	25-30µm	120-150µm
	MPH-XLG-(angle)**	25-30µm	150-180µm
* Angle options include: straight, 15°, 20°, 25°, 30°, 35°, 40°, 45°			
** Angle options include: straight, 15°, 20°, 25°, 30°, 35°			

Micropipet Specifications

- Made from meticulously prepared 1mm O.D. borosilicate glass capillary tubes. Each lot of glass is extensively cleaned and baked before being MEA and endotoxin tested.
- Overall length is approximately 5.4–5.8cm.

Quality Standards

Every pipet is made and checked by a specialized, well trained technician. We offer a 100% guarantee that all pipets meet our standards or we will replace the product.

- CE Marked and FDA Cleared
- ISO 13485 Certified
- Sterilized by gamma radiation
- One-cell MEA and LAL tested

Customization

We know preferences vary; therefore, we will gladly make micropipets to your specifications. For example, you may want a different overall length or a different inner or outer diameter. Simply communicate these preferences when ordering.